


Pipette Management Simplified

Save Time and Money



Optimize your service program

Near- and long-term views of which pipettes are coming due for service or calibration help you better align your lab's workload and deadlines. With PipetteX you can assign different service plans (e.g., calibration vs. preventive maintenance) to each pipette and optimize pipette uptime..

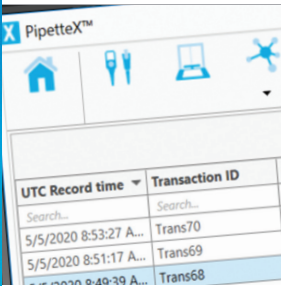

21 CFR Part 11 Compliant

Validate PipetteX for labs that require 21 CFR Part 11 compliance. Managing multiple labs? PipetteX allows you to set policies for only those labs that require them.



Integrate with LIMS and ERP

As your pipette management tool, PipetteX integrates with your LIMS or ERP systems through scheduled exports of calibration data, certificates and routine test reports in structured files.



Implement SmartCheck

Rainin SmartCheck makes periodic pipette performance checks easy. Use PipetteX to schedule routine, between calibration performance checks for specific pipettes or your entire inventory. SmartStand will alert users when pipettes checks are due and SmartCheck will automatically send the results to PipetteX.



PipetteX™ Asset Management

Powerful Software for Busy Labs

In an industry where time = money and time-to-market (or publish) is everything, PipetteX saves time and money. With just a click or two, you know which pipettes need calibration over the next month, quarter or year, and which are out for service. Used in conjunction with Rainin SmartStand™ and SmartCheck, PipetteX simplifies and streamlines pipette management.

PipetteX is flexible, scalable and simple-to-use:

- Automate entering/updating pipette service data for any brand of pipette.
- Improve compliance by automating service and routine test email reminders.
- Automatically capture routine test reports and calibration certificates.

Save Time Managing Pipettes

Optimize Performance and Recordkeeping



Flexible, scalable pipette management

With PipetteX, METTLER TOLEDO has created a simple, yet robust pipette management system that easily adapts to your unique needs.

- Use PipetteX to organize and optimize pipette calibration and recordkeeping – from one lab to multiple labs in many locations.
- Add SmartStands to bring pipette management right to the bench – SmartStand displays each Rainin XLS pipette’s current calibration status, alerts users when a pipette is due for service, and transfers information back and forth between PipetteX and the RFID chip inside every XLS pipette.
- Add SmartCheck to implement routine pipette verification into your workflow – PipetteX can remind users when any pipette is due for performance verification and record the results in the pipette’s record.

Features and System Requirements

Version	Starter*	Site	Unlimited
Catalog number	30475911	30475913	30475915
Pipettes	200	1500	Unlimited
Connected SmartStands	5	200	Unlimited
Connected SmartChecks	5	200	Unlimited
Calibration certificate upload	Manual Upload	Automatic**	Automatic**
Installation type	Single client	Single client/server	Single client/server
Audit trail & e-signature	n/a	Optional with 21 CFR Part 11 license (30587187)	Optional with 21 CFR Part 11 license (30587187)

*Included with SmartStand **Check for local availability

	OS	x32	x64	Edition
Client	Windows 7	Not supported	Not supported	–
	Windows 8.1 with update	Not recommended	Not recommended	Pro, Enterprise
	Windows 10 with update 1709 and above	Conditionally supported	Supported	Pro, Enterprise, Education
Server	Windows Server 2012 R2	–	Not supported	–
	Windows Server 2016	–	Supported	Essentials, Standard, Datacenter
	Windows Server 2019	–	Supported	Essentials, Standard, Datacenter
SQL	SQL Server 2012	Supported	Supported	Enterprise, Standard, Express
	SQL Server 2014	Supported	Supported	Enterprise, Standard, Express
	SQL Server 2016	–	Supported	Enterprise, Standard, Express
	SQL Server 2017	–	Supported	Enterprise, Standard, Express
	SQL Server 2019	–	Supported	Enterprise, Standard, Express

Note: SQL Version below 2012 is not supported by PipetteX – Not Available

Rainin, Pipetting 360®, SmartCheck, PipetteX and XLS are trademarks of Mettler-Toledo Rainin, LLC.

METTLER TOLEDO Group

Local contact: www.mt.com/contacts

www.mt.com/PipetteX

For more information

Subject to technical changes

© 08/2020 METTLER TOLEDO. All rights reserved

17701045 Rev A